

Two-Stage Rotary Vane Pump, Type PK 2 DC

Very small, highly efficient and chemically resistant rotary vane pump for laboratory applications.

Special characteristics:


- free of non-ferrous metals
- media contacting parts made of chemically resistant material
- all elastomers in Viton
- vanes of special plastics
- overload protected motor, IP 54 splash-proof
- water vapor tolerance 33 mbar

Ranges of application:

- vacuum distillation
- vacuum filtration
- vacuum drying
- refrigeration technique etc.

Scope of delivery:

Rotary vane pump ready-to-use with oil charging, centering ring, clamping ring, motor protection, main power switch, mains connection cable and plug.

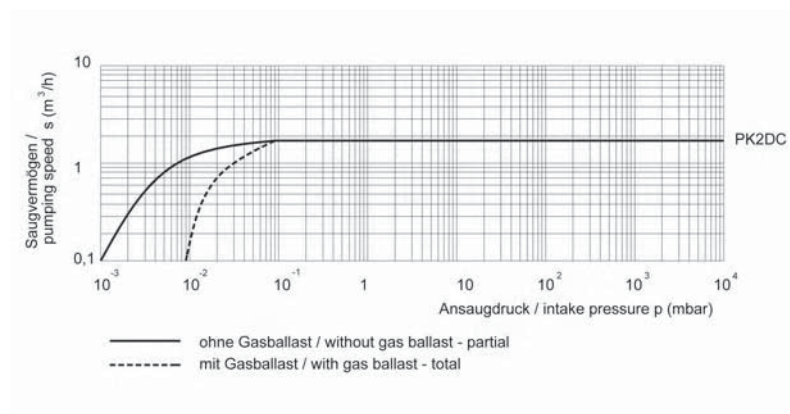

Technical Features

Type	Ult. pressure w/o. gas ballast total mbar	Ult. pressure w/o. gas ballast partial mbar	Ult. pressure w. gas ballast total mbar	Pumping speed pneurop 50/60 Hz m ³ /h	Pumping speed l/min	Dimensions (W/D/H) mm	Weight kg	Motor power W	Oil charge ml	Connections suction/ pressure side
PK 2 DC	1x10 ⁻²	1x10 ⁻³	0.5	1.8/2.2	30/36	300/120/150	8	120	250	DN 16 KF

Ordering Information

Type	Mains supply V / Hz	Mandatory accessories Connection cable	PU pcs.	Order-No.
PK 2 DC	230 / 50/60 CEE	no	1	302010
PK 2 DC	230 / 50/60 UK	no	1	302010-05
PK 2 DC	230 / 50/60 CH	no	1	302010-06
PK 2 DC	115 / 50/60 US	no	1	302010-02

22



Tel.: +49(0)3677 604-0
Fax: +49(0)3677 604-110

E-Mail: info@ilmvac.de
url: www.ilmvac.de



Tel.: +1(0)888 673-0942
Fax: +1(0)757 966-9239

E-Mail: r.askew@ilmvac.com
url: www.ilmvac.us

supplyLAB

www.supplylab.pt
geral@supplylab.pt

Cacém Park - Edifício 9
Estrada de Paço de Arcos nº88
2739-512 Agualva Cacém
T +(351) 21 4278700
F +(351) 21 4278709